


Forced-air Oven H3014

For heat reversion tests and
for general use in laboratories



IptData

- Easy to operate desk-top units
- High accuracy of temperature control
- Temperature range up to 200 °C
- Volumes 56 l or 136 l appropriate for:
 - pipe pieces up to $\varnothing = 200$ mm
 - pipe segments of $\varnothing > 200$ mm

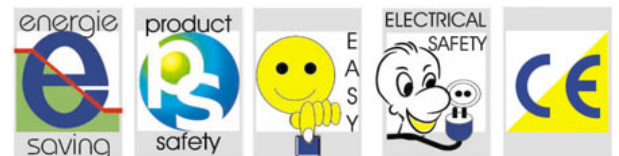

Heat changes dimensions

You expect high temperature constancy

The precise temperature control in combination with a well-contrived air circulation system guarantees a very high spatial and temporal temperature constancy. The circulation system can be controlled in two levels and switches off automatically for approx. 20 sec. after having opened the front door.

Easy handling and safe usage

An integrated time switch serves for an automatic switching off after a preselectable time (up to 999 min.). But as well, the operator can switch the unit on and off manually. All data (temperature, test time) are controlled by a microprocessor. Test time and actual temperature are displayed on a well-arranged control panel.

In your workshop, Ipt is always the right partner for you.

ipt

